
**OPERATION INSTRUCTION
MODE D'EMPLOI
BEDIENUNGSANLEITUNG
INSTRUCCIONES DE EMPLEO
ISTRUZIONI PER L'USO
BRUKSANVISNING**

Swiss LithoClast® Master


Vario ultrasound handpiece
Pièce à main Vario ultrasons
Vario-Ultraschallhandstück
Mango de ultrasonidos Vario
Manipolo agli ultrasuoni Vario
Vario ultraljudhandenhet

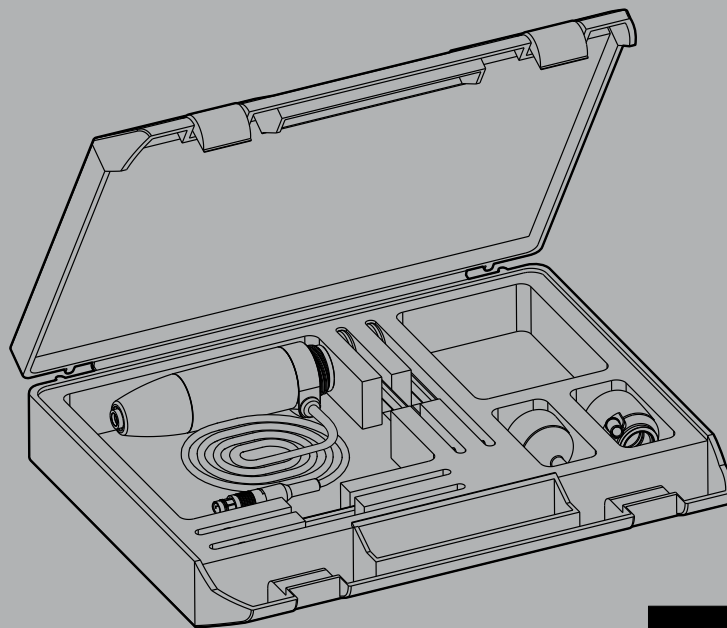ADDITIONAL SAFETY ADVICE FOR USA:

Caution! Federal (USA) law restricts this unit to sale by or on the order of a physician

EMS 
ELECTRO MEDICAL SYSTEMS

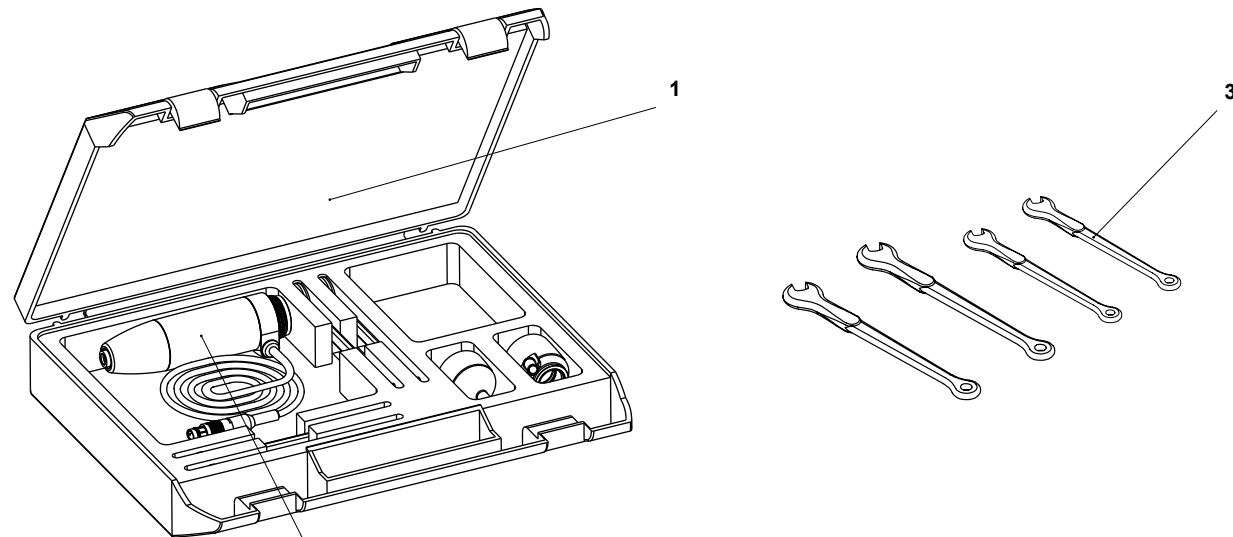COMPONENTS AND CONTROL ELEMENTS

- 1 Ultrasound handpiece Set
- 2 Ultrasound handpiece
- 3 5 mm wrench
- 4 Angled suction connector for combination use ultrasound / pneumatic
- 5 Straight suction connector for ultrasound only use
- 6 Electrical connector
- 7 Protective sterilization cap for the electrical connector
- 8 Protective cap for the Ultrasound Handpiece
- 9 Ultrasound Probe

EMS supplies this product in different configurations.

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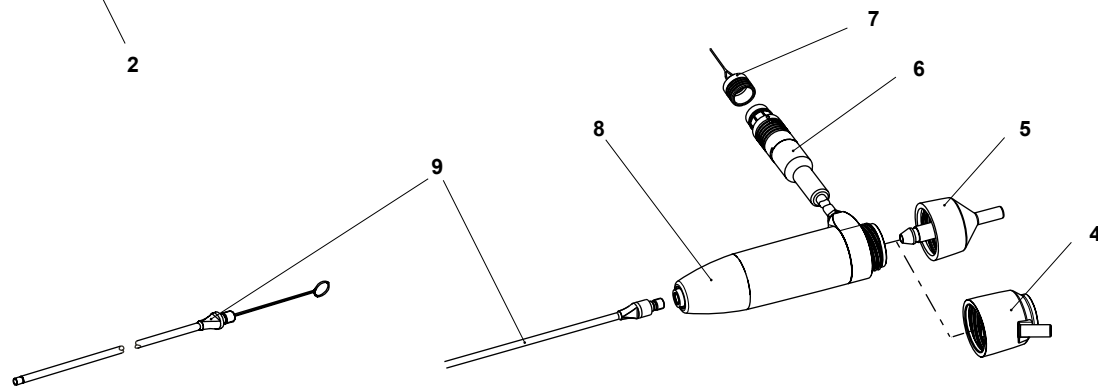
DEAR CUSTOMERS,

Thank you for buying a new EMS product. It meets the highest quality and safety standards.

This Vario ultrasound handpiece is used with the Swiss LithoClast Master for the fragmentation of urinary tract calculi in the kidney, ureter, and bladder.

The Vario handpiece represents an important improvement for performing ultrasound lithotripsy with the Swiss Lithoclast Master. Besides higher effectiveness, improved stability, less weight and smaller size, it provides optimized suction capabilities when working with ultrasound only (not combining ultrasound and pneumatic lithotripsy).

The Vario handpiece is compatible with Lithoclast Master control units which are equipped with software version "rel. 3.0". Control units having a previous version of software installed require an upgrade to enable controlling the Vario handpiece.

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PLEASE READ THIS BEFORE BEGINNING WORK!

About this manual

Note that the English version of this manual is the master from which translations derive. In case of any discrepancy, the binding version is the English text.

These operating instructions are to ensure the correct installation and use of this product.

Please read these operating instructions carefully as they explain all of the most important details and procedures. Please pay special attention to the safety precautions.

Always keep this instruction close at hand.

To prevent injury to people and damage to property, please heed the corresponding directives. They are marked as follows:

-  **Danger**
Risk of injury
-  **Caution**
Risk of property damage or environmental harm
-  **Please note**
Useful additional information and hints
-  **Prohibited**
-  **Authorized**

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We would be pleased to answer your questions or listen to your suggestions. We do, of course, provide support in case of technical problems. Please call directly your approved EMS repair center or your dealer.

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EMS

CONTENTS

Preparing for treatment.....	2
Assembly of the handpiece for ultrasound only use	4
Connecting the ultrasound handpiece to the unit.....	6
Assembly of the ultrasound handpiece for ultrasonic / pneumatic combination use	8
Connecting the ultrasonic / pneumatic handpiece combination to the unit.....	12
Application	14
Coding system for the ultrasound and pneumatic probe sets.....	14
Selection of lithotripsy modality	16
Changing the handpiece assembly from combination mode to ultrasound mode only	22
Function testing	24
Operating the handpiece	24
Cleaning, disinfecting and sterilizing	28
Preparation of the ultrasound handpiece for reprocessing	28
Preparation of the combination handpiece for reprocessing	30
Maintenance.....	34
O-ring	34
Safety precautions	36
Storing of the product	40
Product disposal	42
EMS-Service.....	44
Symbols	46
Technical data	48
Troubleshooting	52

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PREPARING FOR TREATMENT

 Clean, disinfect and sterilize all the components before first use and after as described in the "Reprocessing Instructions" manual. Assemble the application parts under sterile conditions.

Before treatment, perform a function test according to the instructions given in the respective paragraph.

 Sterilized backup probes should be available in case of a probe failure during treatment. Should a probe break distally, sterile grasping forceps should be used to remove probe pieces from the urinary tract.

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

**Assembly of the handpiece for ultra-
sound only use**

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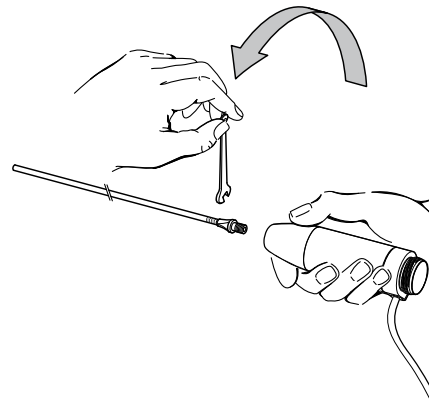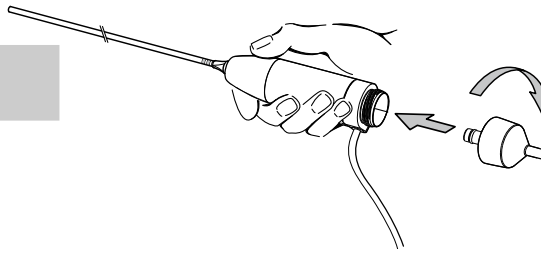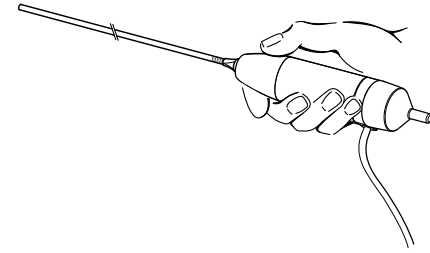
- ⚠ The ultrasound probe must be tightened onto the ultrasound handpiece using the supplied wrench. Keep the wrench sterile. If the probe loosens during treatment, retighten the probe with the wrench. Even a slightly loose probe connection will cause a total loss of fragmentation performance.
- ⚠ To thread the suction connector in properly, align it carefully with the handpiece and apply a light axial pressure while screwing it in place.

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Connecting the ultrasound handpiece to the unit

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
Ansluta ultraljudhandenheten till enheten

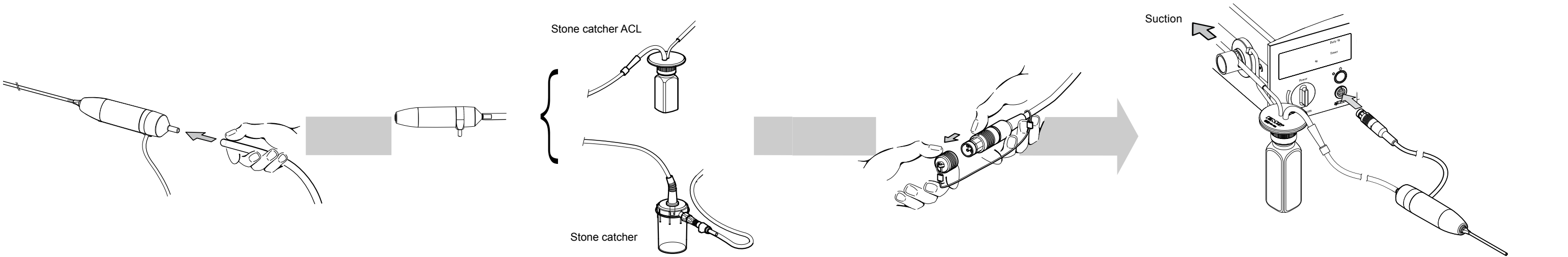For further details regarding the preparation of the Stone catcher ACL, please refer to the operating instructions of the Stone catcher ACL.

Connect the long tube of the sterile single use Stone catcher to the suction connector of the ultrasonic handpiece.

The electrical connectors must be perfectly dry before connecting them to the unit.

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

Assembly of the ultrasound handpiece for ultrasonic / pneumatic combination use

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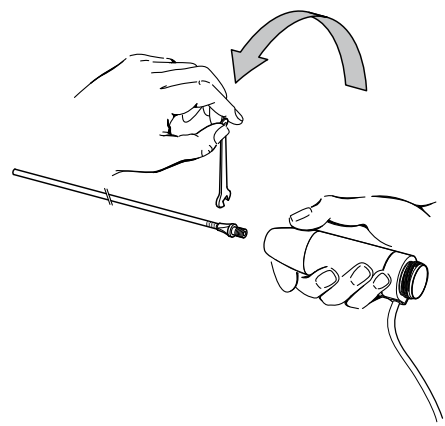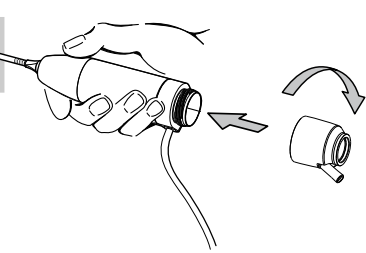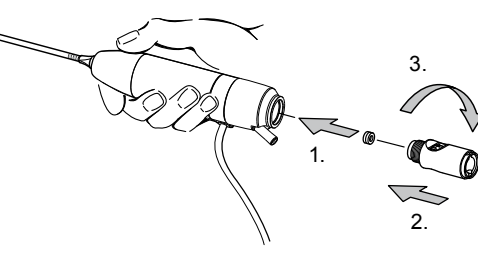
⚠ The ultrasound probe must be tightened onto the ultrasound handpiece by using the supplied wrench. Store the wrench sterile. If the probe loosens during treatment, retighten the probe by using the wrench. Even a slightly loose probe connection will cause a total loss of fragmentation performance.

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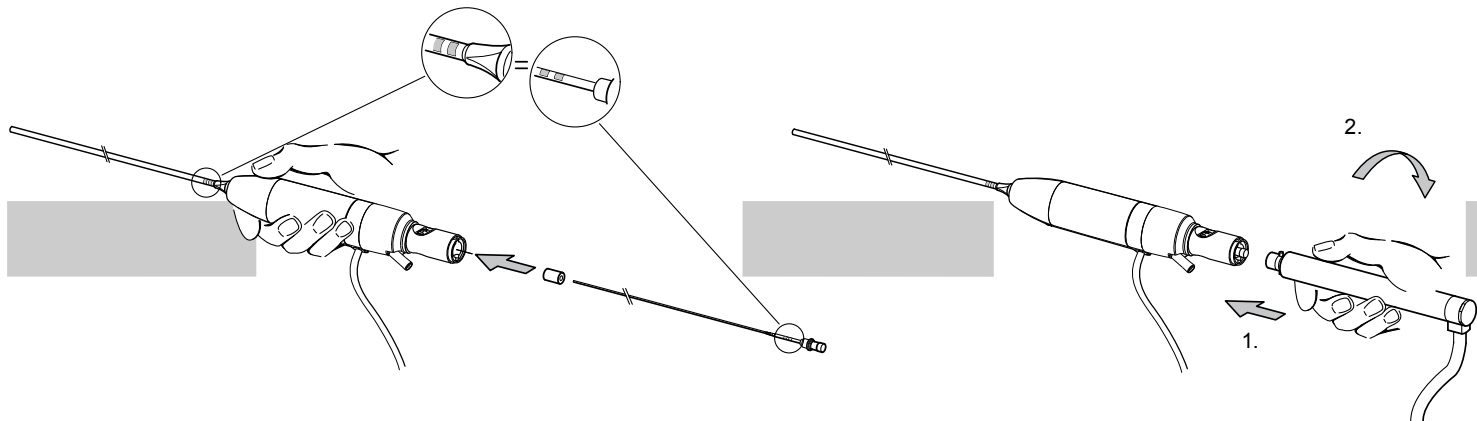
Adjust the length of the pneumatic probe by turning the wheel of the adjustment interface so that the tip of the pneumatic probe is flush with the tip of the ultrasound probe.

 If the pneumatic probe cannot be adjusted flush with the ultrasound probe, this indicates the wrong ultrasound / pneumatic probe combination is being used. Combination probes are marked by an identical ring pattern on the proximal shaft of both probes.

 Both lithotripters work only in direct contact with calculi. The length adjustment of the pneumatic probe relative to the ultrasound probe is important, when using both lithotripters simultaneously.

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


UTILISATION

ANWENDUNG

Coding system for the ultrasound and
pneumatic probe sets

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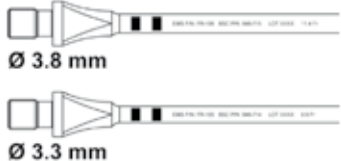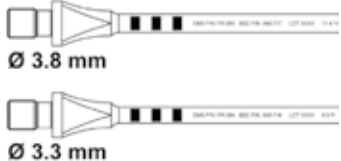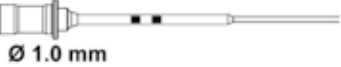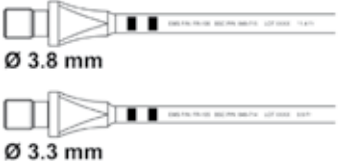
Selection of lithotripsy modality

The following table shows which ultrasound probe size and pneumatic combination probe size is used for various stone types and locations:

	Lithotripsy Mode	
Stone Type	Ultrasound only	Pneumatic & Ultrasound combination
Hard kidney stones through percutaneous access (monohydrate, cystine stones)	Ultrasound probe: 3.3 mm or 3.8 mm Use ultrasonic lithotripsy only for hard stones after initial fragmentation with the pneumatic / ultrasonic combination (see directions for intraoperative change between combination handpiece setup and ultrasound only)	Ultrasound probe: 3.3 mm or 3.8 mm Pneumatic probe: 1.0 mm Pneumatic probe 1.3 mm (optional) for use with 3.8 mm ultrasound probe only
Soft kidney stones (struvite, matrix stones)	Ultrasound probe: 3.3 mm or 3.8 mm	
Ureteral stones	Ultrasound probe: 1.5 mm	
Hard bladder stones (monohydrate stones)	Ultrasound probe: 3.3 mm or 3.8 mm This method can be used for fragment removal following the pneumatic fragmentation	Ultrasound probe: 3.3 mm or 3.8 mm Pneumatic probe: 1.0 mm Pneumatic probe 1.3 mm (optional) for use with 3.8 mm ultrasound probe only
Soft bladder stones	Ultrasound probe: 3.3 mm or 3.8 mm	
Small kidney stones	Ultrasound probe 1.9 mm for use with mininephroscopes with 6 Fr. work channel size	

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


**Changing the handpiece assembly
from combination mode to ultrasound
mode only**

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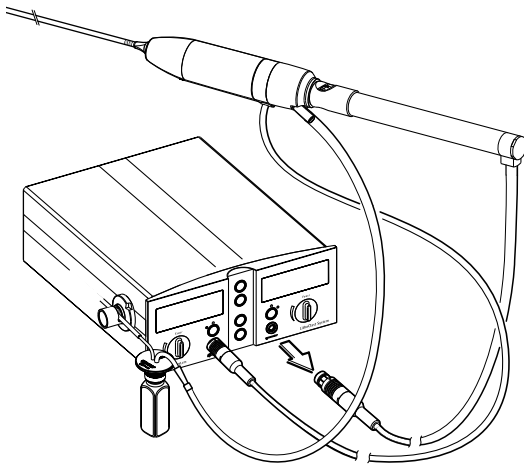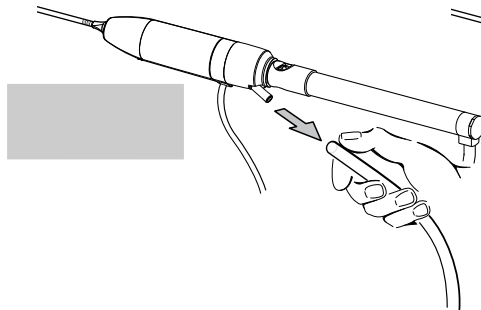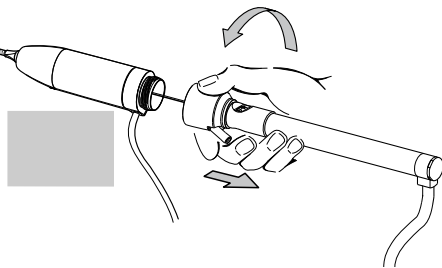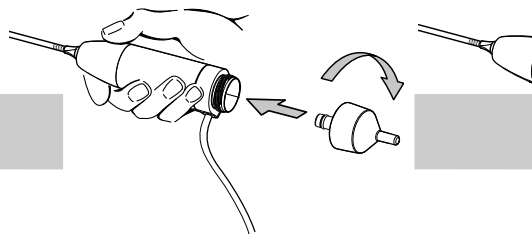
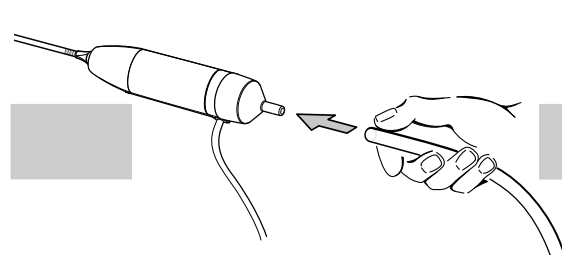
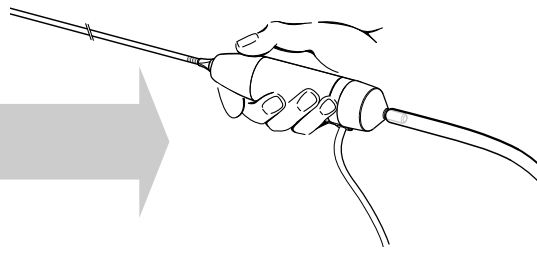
👉 The pneumatic connector stays connected if the pneumatic handpiece with probe is removed temporarily.

👉 To change from combined ultrasonic / pneumatic lithotripsy to stand alone ultrasound lithotripsy, it is necessary to replace the angled suction connector by the straight one.

⚠ Always use the straight suction connector when working with ultrasound only.

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Function testing

 Prior to using the handpieces, a function test should be performed. Check that aspiration is always connected and operational.

 Immerse the tip of the ultrasound probe into a cup with sterile water. Activate the ultrasound and make sure that suction is functioning and that ultrasound vibration is present (presence of a continuous noise) at the probe tip.

Operating the handpiece

 Do not touch the ultrasound probe when the ultrasound is activated. The ultrasound vibrations in the probes may cause instant burns by frictional heating.

 Danger of urothelial and mucosal perforation! There is a risk of injuring the mucosa of the kidney and bladder, which lies behind the stone. Always work under endoscopic visualization and keep the pneumatic probe tip flush with the ultrasound probe tip or the LithoVac probe when using the combined modes.

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 Fragments can block suction through the ultrasound probe. This will lead to a loss of the ultrasound suction capabilities and to heat-up of ultrasound probe and handpiece. Check at regular intervals that the suction is functional.

 The ultrasound handpiece and probe error detection is only possible when working with 100% "Duty cycle".

 For further details on how to use this product in combination with the Swiss Lithoclast Master, please refer to the instruction manual of the Swiss Lithoclast Master.

 Check that the suction connector is tight on the handpiece to maintain a good seal on the suction system.

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 Please conform to the recommendations of the "Reprocessing Instructions" manual delivered with your product regarding procedure of cleaning, disinfecting, packaging and sterilizing of the components.

 The handpiece must be fully disassembled for each sterilization.

 Make sure the probe is replaced every 5 sterilization cycles.

 Please follow the applicable present day sterilization regulations in your country.

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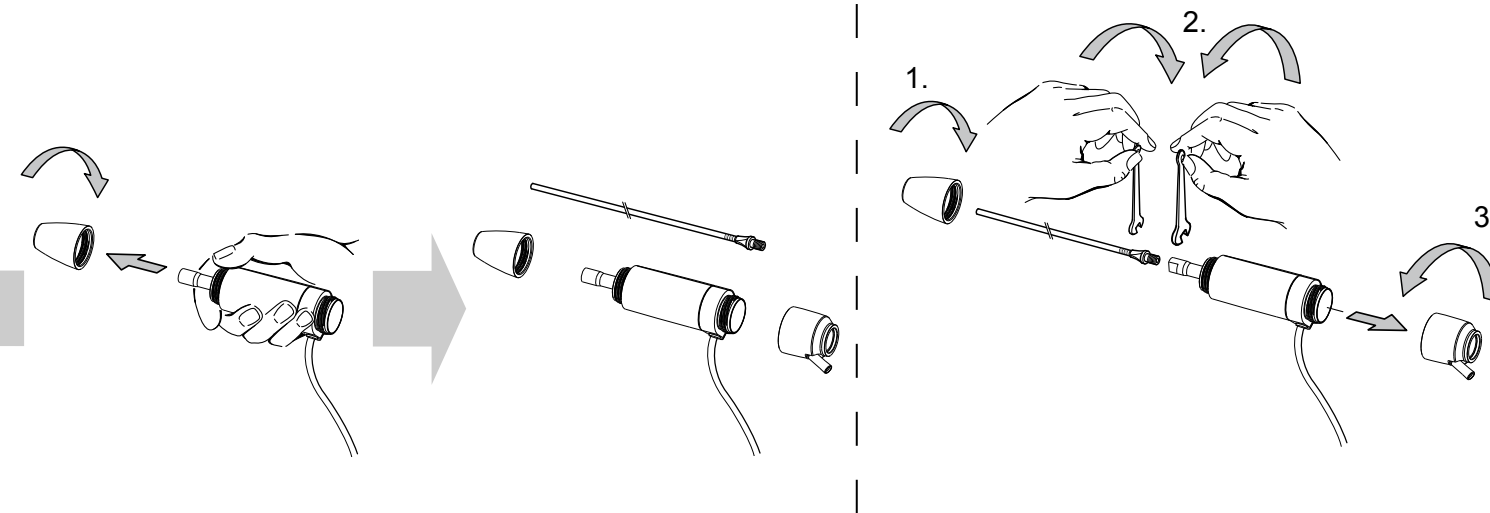
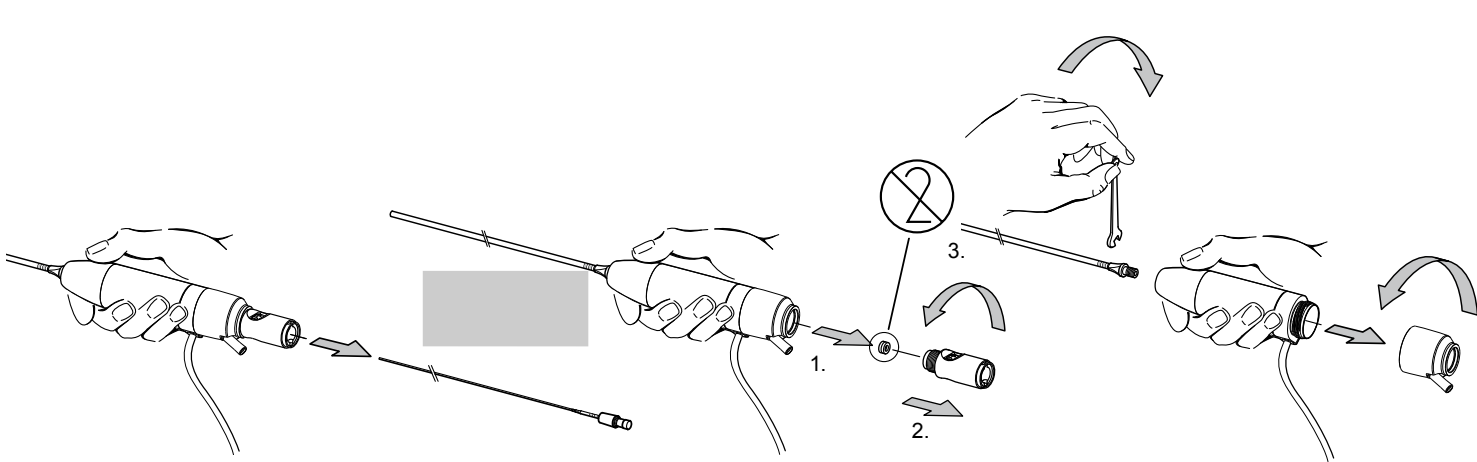
⚠ In each sterilization procedure, please be sure that the silicone seal has been replaced and discarded.

⚠ Check all wear parts, such as probes, probe guides, etc. regularly, and replace if necessary.

⚠ The ultrasound handpiece must be fully disassembled for each sterilization.

⚠ Ensure sterilization cap is connected to each hose before processing.

👉 If the ultrasound probe is hard to untighten from the handpiece, use the 5 mm wrench and the 8 mm wrench as shown to unlock.

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SAFETY PRECAUTIONS

EMS and the distributor of this product accept no liability for direct or consequential injury or damage resulting from improper use, arising in particular through non-observance of the operating instructions, or improper preparation and maintenance.

 Use for the intended purpose only: Before using the product, make sure that you have studied the operating instructions. This also applies to any equipment used with this product. Failure to observe the operating instructions may result in the patient or user suffering serious injury or the product being damaged, possibly beyond repair.

 Risk of explosion: Do not use this product in the presence of inflammable anaesthetics or gases.

 Use of the device is contra-indicated in patients with the following:

- Active bleeding disorders
- Solitary functioning kidney
- Creatinine greater than or equal to 3 µg %
- During pregnancy
- Stricture and obstruction problems
- With an electrical stimulator implanted e.g. pacemaker.

 This product must be used only by trained and qualified personnel.

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Before each use, always clean, disinfect and sterilize the different pieces and accessories of this product. Please refer to the information provided in the "Reprocessing Instructions" manual delivered with the product. Non-sterile pieces and accessories may cause bacterial or viral infections.

Always examine the product for damage before commencing treatment. Damaged accessories or a damaged unit must not be used and must be replaced. Use original EMS spare parts and accessories only.

Have this product exclusively repaired by an approved EMS repair center.

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STORING OF THE PRODUCT

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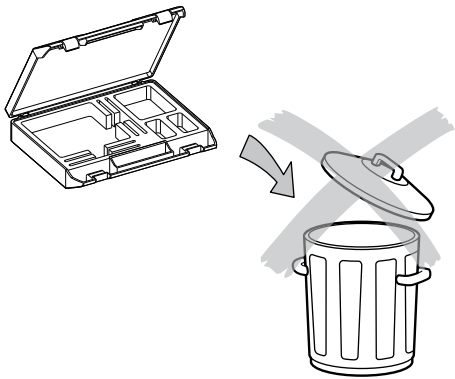
 Keep the original packaging until the product is to be disposed of permanently. It can be used for shipping or storing at any time.

Should you wish to store the product for an extended period of time:

- Proceed as described in the chapter "Cleaning, disinfecting and sterilizing",
- Pack the product and all accessories in the original packaging.

Please refer to the storage and transport conditions in the chapter "Technical data".

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PRODUCT DISPOSAL

The product must not be discarded in domestic household waste.

Should you wish to definitively dispose of the product, please comply with the regulations which apply in your country.

WARRANTY

The warranty is valid from the date of purchase of your product for the period indicated in the "Registration Card".

Damages due to non-adherence to the operating instructions or wearout of parts are excluded from warranty.

 The warranty of your product will be cancelled if you try to open it.

ACCESSORIES

Accessories are available from EMS or any authorized dealers. Please contact your customer service directly.

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EMS-SERVICE

Should your product need additional servicing or repairs, please send it to your dealer or our approved EMS repair center.

In the case of non-authorized repairs or damaged due to non-adherence to the operating instructions, EMS accepts no liability whatsoever. This will also void the warranty.

It is best to ship your product in the original packaging. It protects your product against damage during shipment.

 Before dispatching your product including all accessories, please clean, disinfect and sterilize as described in the operating instructions.

 Risk of transmitting micro-organisms ! To protect the personnel of the repair center and for safety reasons during transport and shipment, all products and accessories returned to the factory for repair or revision must be cleaned, disinfected, rinsed and sterilized in accordance with the instructions manual. Repair can be refused for product or accessories received in a contaminated condition.

When sending your product directly to the approved EMS repair center, please include the name of your dealer. This simplifies processing for us.

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TECHNICAL DATA

DESCRIPTION	
Manufacturer	EMS SA, CH-1260 Nyon, Switzerland
Model	Vario ultrasonic handpiece for Swiss Litho-Clast® Master
Ultrasound frequency	23.2 – 26.4 kHz
Classification EN-60601-1	Class I Applied part BF IP67 Mode : Continuous operation
Classification MDD 93/42 EEC	Class IIb
Weight	270 g
Dimensions	Length : 132-138 mm Width : 31.5 mm
Operating conditions	+10°C to +40°C 30% to 75% relative humidity
Validated life time	100 treatments
Storage and transport conditions	-10°C to +40°C 10% to 95% relative humidity 500 hPa to 1060 hPa air pressure

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TROUBLESHOOTING

 For all codes displayed on the unit may, please read the instruction manual of the unit.

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EMS Corporation
EMS Medical GmbH